

Welcome to the New Age of Digital X-Ray

A complete digital x-ray solution designed for low volume practices from the world's leader in digital imaging.

The FCR Prima™ is the newest member of Fujifilm's digital x-ray family of products customized for private practices. It has been designed to meet the needs of low imaging volume practices without compromising image quality. FCR Prima is one of the most compact and lightweight CR systems available. Taking up less than 2.5 sq. ft, FCR Prima can be installed virtually anywhere.



Exceptional image quality, every time.


FCR Prima scans at 10 pixels/mm (5 LP/MM) giving you sharp image detail. Only Fujifilm offers Image Intelligence™ – the culmination of over seventy years of image processing experience for stable, optimized, and high-quality images every time.



FCR Prima Specifications:

Standard Components	FCR Prima Image Reader (Model: CR-IR 391RU) AC power cord
Other System Components	FCRView Image Management Workstation, FCRView-Prima Cart (optional), Dry Imager: DRYPIX 2000, 4000, 7000 (optional)
Supplies	Imaging Plate ST-VI: 14" x 17", 14" x 14", 10" x 12", 8" x 10", 24 x 30 cm, 18 x 24 cm, 15 x 30 cm IP Cassette Type CC: 14" x 17", 14" x 14", 10" x 12", 8" x 10", 24 x 30 cm, 18 x 24 cm, 15 x 30 cm
Processing Capacity	Up to 29 Imaging Plates/hr.
Time to Image Display	Min. 23 sec.
Dimensions (W x D x H)	23.6" x 15.7" x 30.7"
Weight	154lbs.
Power Supply Conditions	Single phase 50-60Hz AC100-240V ±10% 5A (max)
Environmental Conditions	Operating Conditions: • Temperature: 59-86° F • Humidity: 15-80%RH (No dew condensation) • Atmospheric pressure: 750-1060hPa

Dimensions:

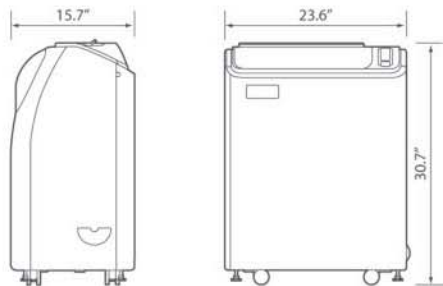

For more information, contact Fujifilm at 1-866-879-0006 or www.fujiprivatepractice.com



Authorized Fujifilm Dealer:

Essential Enterprise Solutions, Inc.
3686 Old Keystone Road
Tarpon Springs, FL 34688
(727) 938-3374 - Office
(727) 499-7814 - Fax
<http://www.eesimed.com>

Flexible reading in multiple cassette sizes.

The FCR Prima cassette supports multiple sizes for a wide range of studies, including chest, lumbar spine, panoramic (15 x 30cm), and extremities.

Imaging Plate ST-VI

14" x 17", 14" x 14",
10" x 12", 8" x 10",
24 x 30 cm, 18 x 24 cm,
15 x 30 cm



IP Cassette Type CC

14" x 17", 14" x 14",
10" x 12", 8" x 10",
24 x 30 cm, 18 x 24 cm,
15 x 30 cm



Durability you can rely on.

FCR systems have proven reliability to meet the rigors of the private practice, thanks to the intelligent design and built-in redundancy features that prevents problems before they can occur. FCR's minimal contact cassette, dust cover and internal optical scanning brush prevents dust artifacts and maximizes uptime.



Shown with optional cart

FUJIFILM